




Traction Models

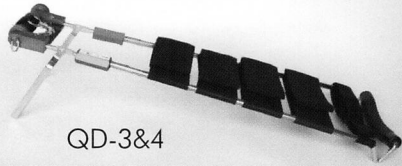

QD-3&4

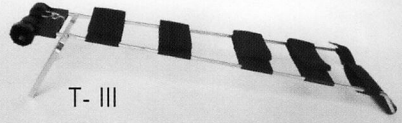

T- III

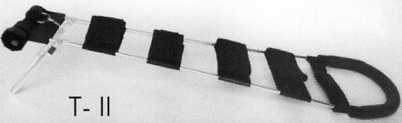

T- II



QD Leg Traction Kits

The QD line of splints kits has several unique and proprietary features. By centering the Hex Lock the splint can quickly be adjusted to any size. Once the length has been determined, the intuitive locking mechanism is engaged by turning the Hex Lock to either side. The Hex Lock system has been designed to work in the worst of weather conditions. No concern with thread direction and no problems with dirt or corrosion which cause other splints to malfunction.

The adult splint will adjust to accommodate from pediatric to adult size range.

All splints come as a complete kit with straps, ankle hitch and carrying case included.

All splints are available with either our anodized aluminum ratchet or our injection molded unit.

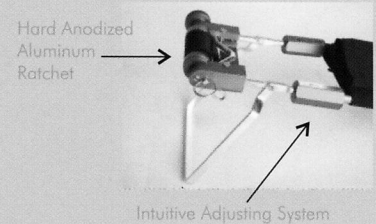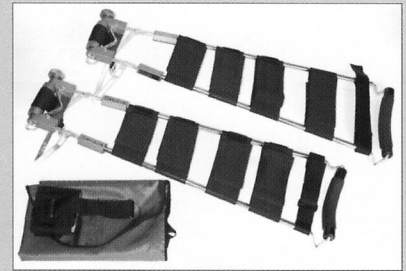

QD 3 & 4 Traction Kits - Adjustable Length

Part Number	Description	Ratchet type
1126524	QD-4 Adult Splint	Aluminum
1126522	QD-3 Pediatric Splint	Aluminum
1126514	QD-4 Adult Splint	Molded
1126512	QD-3 Pediatric Splint	Molded
1126726	QD 3&4 Combination	Aluminum
1126716	QD 3&4 Combination	Molded

TIII Traction Leg Splint kits (Hare type) Fixed Length

Part Number	Description	Ratchet type
1123526	TIII Splint Kit Adult	Aluminum
1123524	TIII Splint Kit Pediatric	Aluminum
1123516	TIII Splint Kit Adult	Molded
1123512	TIII Splint Kit Pediatric	Molded

TII Traction Leg Splint (Thomas Type)

Part Number	Description	Ratchet type
1123616	TII Thomas Type Kit Adult	Molded
1123612	TII Thomas Type Kit Pedi	Molded

Thomas Splint- Frame Only

Part Number	Description	Ratchet type
1223816	Frame only Adult 52" LOA	None
1223812	Frame only Child 31" LOA	None

Specs	Length Extended	Length Closed
QD-4	54"	32"
QD-3	45"	29"

Replacement Ratchets - with straps	
1124500	Aluminum
1124300	Injection Molding

Replacement Leg Slings- fits all models	
1123916	Adult
1123912	Child

Ankle Hitch - Universal	
1124200	Fits all Sizes

Replacement Carrying Cases	
1125010	QD-4 Adult Case
1125008	QD-3 Pediatric Case
1125026	TII and TIII Adult Case
1125022	TII and TIII Pediatric Case

Wire Ladder Splints


This staple of the EMS Industry is alive and doing well. Still manufactured in large numbers, the Faretec Wire Ladder Splint is coated with a thick epoxy coating to cover all edges. They are pliable and can be formed to any shape.

This is one of the most economical forms of utility splint on the market today.

They are packaged and sold by the dozen.

Dimensions: 31" x 3.75"

Part # 1211300

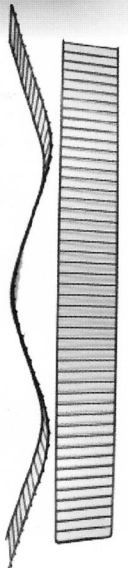

Faretec RediSplint

New in 2002, the Faretec RediSplint is the latest splint in our product line. The RediSplint is made from an alloy that is stiffer and more durable than other utility splints now on the market. With thicker padding on both sides, the RediSplint has been designed for increased comfort. In addition, the padding is non-absorbant, which may allow the splint to be cleaned and reused. The splint folds up into a compact size for economical use of storage space. Currently the RediSplint is packaged and sold in a box of four. Other sized RediSplints will be available soon.

Dimension : 36" x 4"

Part # 1211302

